


# ENOCHS REGENCY 450

## OB/GYN Examination Table

**PATIENT LOAD: 600 LBS.**

evenly distributed on the table's main couch

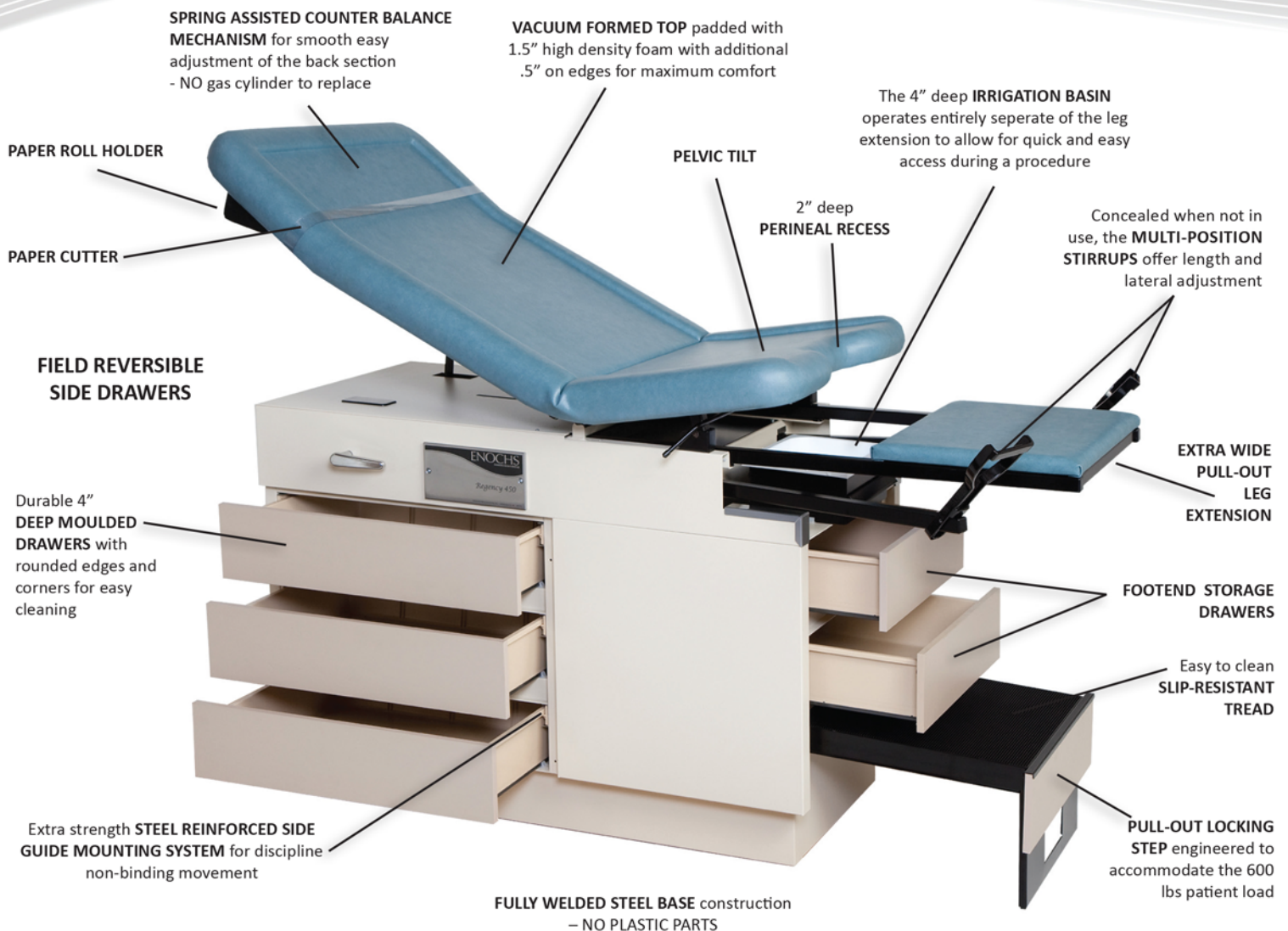

### TABLE FEATURES:

- 2-YEAR Warranty
- Fully Welded all Steel Base and Top Construction - NO PLASTIC
- Spring Assisted counter balance allows smooth adjustment of couch back
- Easy to clean 100% Seamless Vacuum Formed Top
- Pelvic Tilt with a 7° and 11° adjustment
- 4" deep Irrigation Basin engineered to move independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers
- Certified by UL to Standard 60601-1 CAN/CSA C22.3 No. 601-1
- ENOCHS Regency 450 replaces the Regency 60-01



The working end of the table features the independent operation between the stirrups, leg extension and irrigation basin. This efficient design gives practitioner fast and easy transition from one to another. The padded leg extension pulls smoothly out from the table to support patient's legs and can be pushed under the seat section, allowing full access to the irrigation basin and perineal recess.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The spring assisted counter balance lift mechanism (NO gas cylinder) allows for quick and convenient adjustment of the couch back to an infinite number of positions from horizontal to a full 80° chair position. The Vacuum formed top, designed for strength and durability has an all steel back and frame (NO wooden or plastic parts). The durable Naugaform® covered top is easy to maintain and provides a long, trouble-free service life.

# Specifications for the ENOCHS REGENCY 450 OB/GYN Examination Table

<b>Warranty:</b>	2-YEAR
<b>Base Construction:</b>	Fully welded all steel base
<b>Back Lift Mechanism:</b>	Horizontal to 80° by spring assisted counter balance
<b>Couch Top:</b>	100% Seamless vacuum formed style All steel frame construction
<b>Perineal Recess:</b>	Moulded, no seams, 10" W x 2" D at center point
<b>Pelvic Tilt:</b>	7°, 11°
<b>Leg Extension:</b>	Extends out 20" from end of seat section
<b>Irrigation Basin:</b>	4" Deep basin
<b>Electrical Outlets:</b>	Not available on this model
<b>Instrument Warmer: with ON/OFF Switch, Circuit Breaker &amp; Reset</b>	Not available on this model
<b>Stirrups:</b>	Multi-positional, lateral and length adjustment
<b>Storage:</b>	<b>3 Side Drawers:</b> 21.5" W x 15" D x 4" H <b>2 Footend Drawers:</b> 15.5" W x 15" D x 4" H <b>Total Storage Space:</b> 3.32 cubic feet
<b>Foot Step:</b>	Extra long surface engineered for 600 lb. patient load
<b>Paper Roll Holder:</b>	Accommodates 18" or 21" paper widths
<b>Accessories:</b>	knee crutches, arm board, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, headend accessory rails, footend accessory rails, caster base
<b>Mains Power Entry: with Circuit Breaker</b>	Not available on this model
<b>Duty Cycle:</b>	<b>Receptacles:</b> n/a <b>Instrument Warming Drawer:</b> n/a
<b>Ratings:</b>	<b>Table:</b> n/a <b>Receptacles:</b> n/a <b>Total Receptacle Load:</b> n/a
<b>Table Dimensions:</b>	<b>Width:</b> 27", <b>Height:</b> 34", <b>Length:</b> 52" (extends to 72")
<b>Table Weight:</b>	268 lbs.
<b>Patient Load:</b>	600 lbs. evenly distributed on main couch
<b>Certifications:</b>	By UL to Standard 60601-1, CAN/CSA C22.3 No. 601-1